

Item Name : Free Fall Apparatus

Item Code : 3067.2008.01


Item Description :

- Micro-magnet holds ball in start position
- Three contact pins ensure start position can be reproduced and act as contacts of a switch that opens when the ball is released, triggering the timing measurement
- When the ball strikes the contact plate, timer stops
- Drop height can be adjusted to a fraction of a millimeter
- Measurement read off a scale on the column

Specifications

- Height scale: 20 to 960mm • Scale divisions: 10mm
- Scale precision: 0.2mm • Balls: Steel, 16 mm dia.
- Dimensions: ~200 x 130 x 1000 mm³ • Weight: ~1.6kg

Item Image :

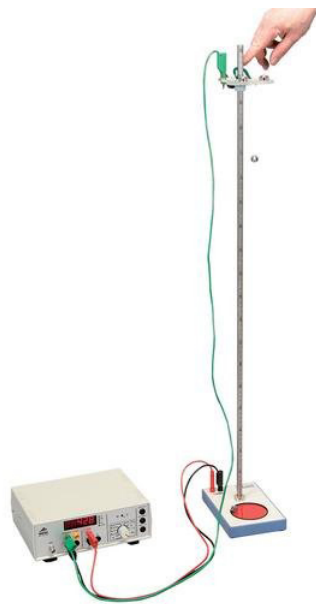

Cat. No.

3067.2008.01